

BeHear[®] PROXY

Wireless Neck Speaker and
Personalized Hearing Amplifier



Surround Sound on Your Shoulders


- Keep ears free with this comfortable, flexible, wireless neck speaker
- Crisp TV dialogue and rich stereo sound just inches from the ears
- Loud and clear wireless mobile calls
- Personalized amplification of live conversations and other world sounds
- Extractable earphones for privacy, or in consideration of others
- Supports two Bluetooth devices simultaneously (e.g., smartphone and audio transmitter)



Suitable for use with or without hearing aids


TV Watching: Exceptional Audio

- Top quality full stereo sound with Alango's SonicVibrance™
- Minimal speaker-to-ear distance eliminates room reverberation, reduces the required sound volume and listening effort
- Easily position this lightweight speaker to rest comfortably on your shoulders
- Can be used with hearing aids

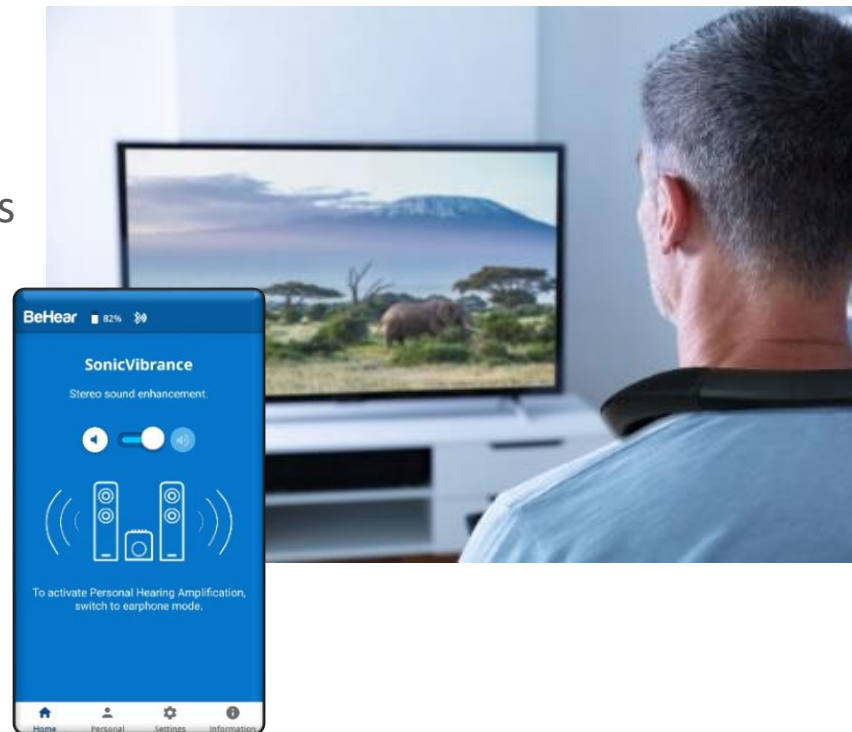

Mobile Calls: Loud and Clear


- Wireless phone calls via Bluetooth
- Call control (accept/reject/end, volume) directly from the neck speaker
- Improves intelligibility with personalized sound, and by slowing down speech of “fast talkers” with EasyListen™
- Extractable earbuds keep conversation private, when desired

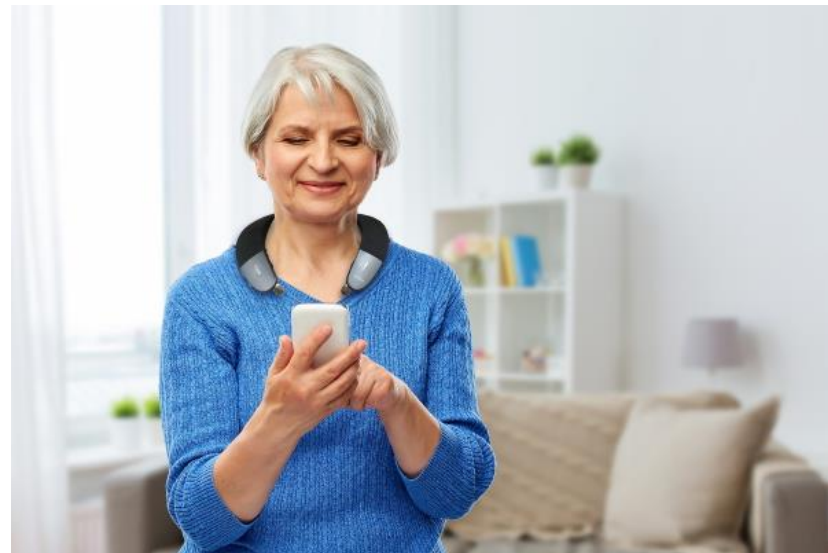

Easy Sound Personalization


1. Pair the BeHear PROXY with a smartphone using the W&H BeHear app.
2. Using the app, perform the simple hearing assessment.
3. Results of the assessment are automatically programmed into the neck speaker.

From now on, all important sounds are amplified and clarified to match the wearer's unique hearing profile.

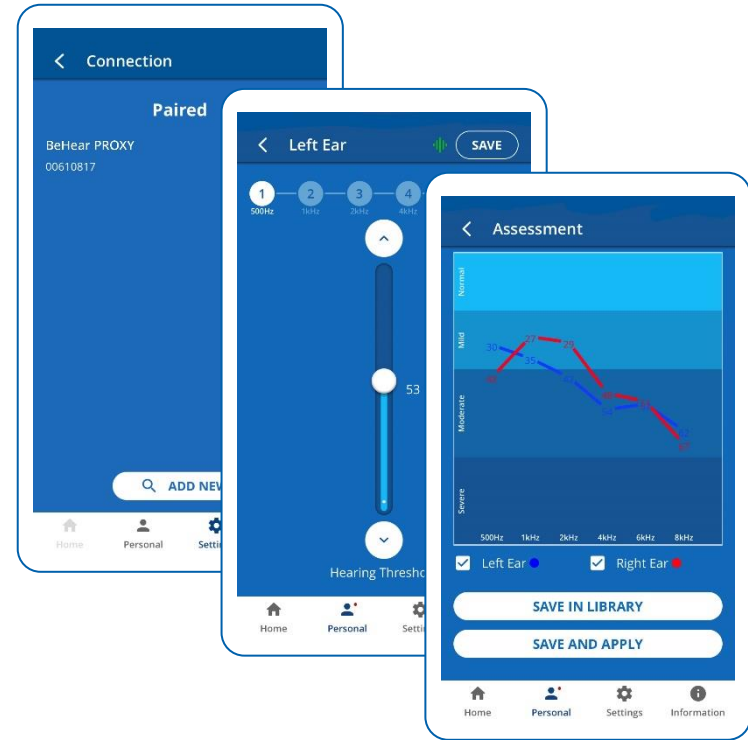

Personalized Hearing Amplifier

- Retractable earbuds convert BeHear PROXY into a headset for private listening and personal hearing amplification
- Speakers are automatically disabled in earphone mode
- Great for private calls, or in consideration of others
- Can be used as a back-up hearing aid

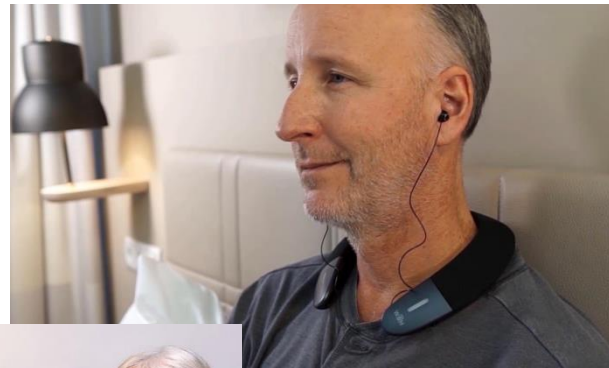


Optimal TV Watching Experience

- Pairs with HearLink® PLUS* Bluetooth audio transmitter
- Long range (up to ~60m / 200ft) enables freedom of movement
- Fully synchronized audio-video: aptX™ low-latency
- Support for both optical and analog audio outputs
- TV audio controls are located on the neck speaker
- For mobile calls: auto pause/resume of audio stream



** Sold separately*

Neck Speaker Functions

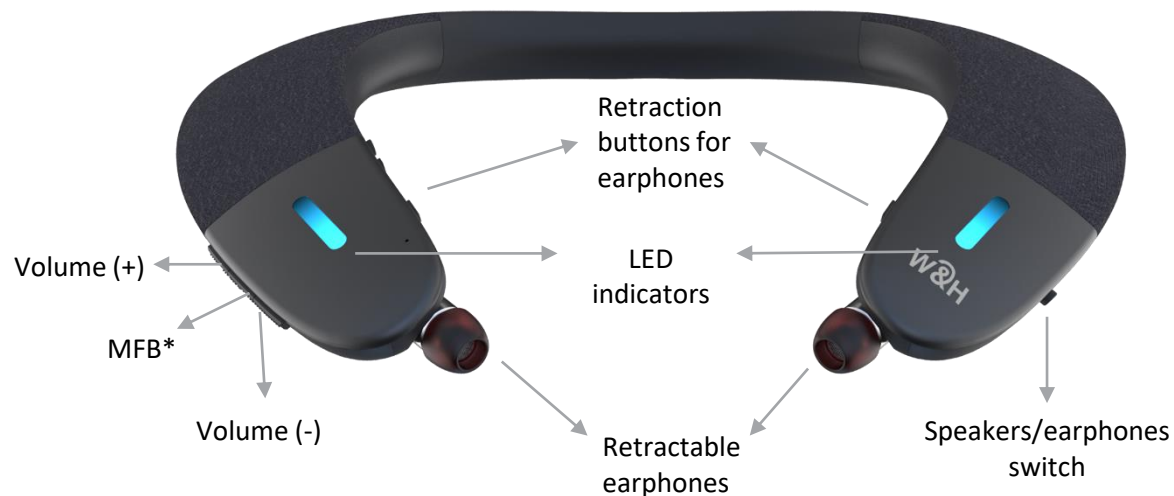

* MFB – Multi-Function Button (pairing mode, answer/reject call)

Technical Specifications

- **Microphone**
 - Type Omni-directional MEMS
 - Frequency range 50Hz – 10kHz
- **Earphones**
 - Frequency range 20Hz - 20kHz
 - Max SPL 110dB
- **Speakers**
 - Frequency range Up to 20kHz
 - Max SPL 110dB
- **Bluetooth**
 - v.4.2
 - aptX Low Latency
- **Battery**
 - Capacity 1050mAh
- **Operation Time**
 - Personal sound amplification up to 50 hours
 - Calls using earphones up to 45 hours
 - Media play (earphones) up to 40 hours
 - Media play (loudspeakers) up to 10 hours
- **Weight**
 - 256g (9 oz.)

Package Contents

Mushroom-Shaped Tips

4x Pairs



Cone-Shaped Tips

3x Pairs



Open Fit Tips

1x Pair



Open Ear Tips

1x Pair



BeHear®
PROXY



Micro USB
Charging Cable



Quick Start
Guide



Thank You

From experts with love:
Sound Enhancement by Alango Technologies, Ltd.

info@wearandhear.com

www.wearandhear.com